

# Aesculap Instruments for tunnel technique

To cover periodontal recessions



Aesculap Dental

# Instruments for tunnel technique




Rounded plate-shaped working tips allow more movement and significantly increase the radius of mobilisation in frontal and lateral sub-periosteal preparation.

The sharpened edge of the instrument tips also allows sulcular as well as tunnel technique preparation (as with the envelope technique).



The design and angle of the working tips significantly facilitates preparation in convex areas.

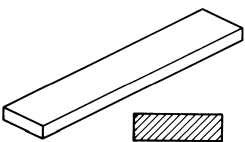


DB095

Original Arkansas grinding stone for grinding on tunnelling instruments  
75 x 25 x 7 mm



JG598

Sterilit® Aesculap oil for grinding stones  
Contents: 50 ccm

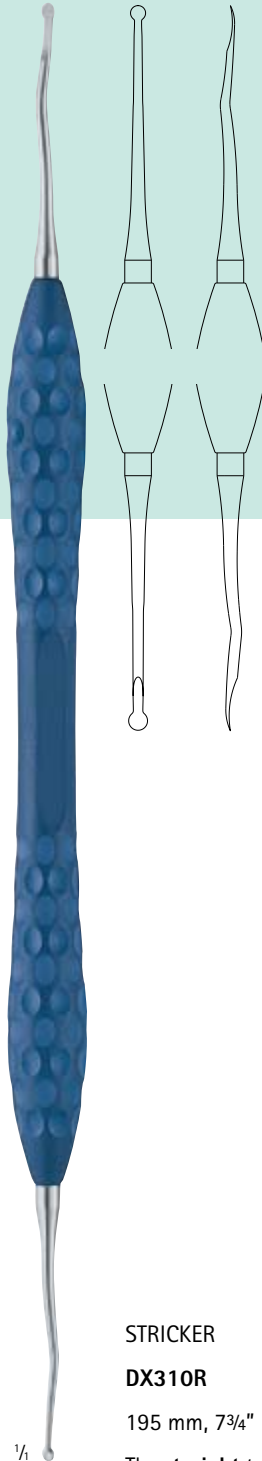


STRICKER  
DX310R

195 mm, 7<sup>3</sup>/<sub>4</sub>"

The **straight** tunnelling instrument is mainly used in the upper jaw, in the frontal and lateral region. This allows large areas of the gingiva to be undermined using both tunnel and sulcular techniques.



# for periodontal recessions

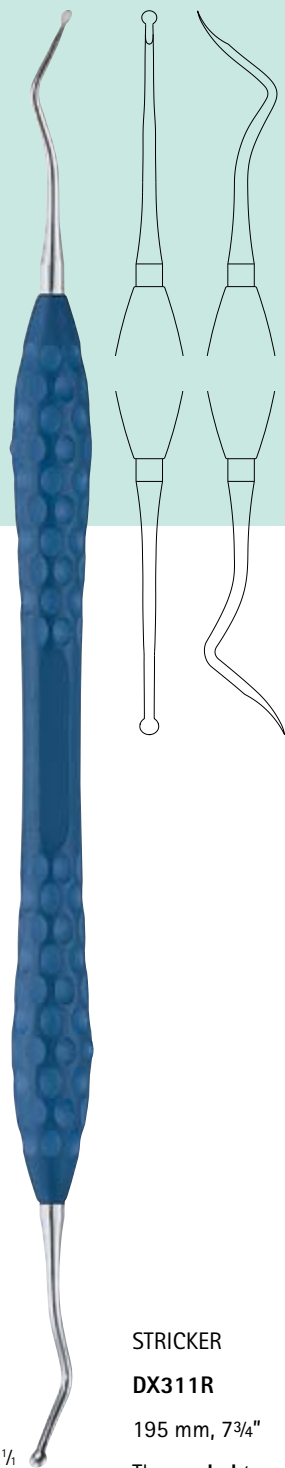

**STRICKER  
DX311R**

195 mm, 7<sup>3</sup>/<sub>4</sub>"

The **angled** tunnelling instrument is primarily used for preparing lateral, buccal regions. It can also be used distally, without tearing the mucosal tissue and also for the lower jaw.



**BB045T**

145 mm, 5<sup>3</sup>/<sub>4</sub>"

Titanium scalpel handle for micro-surgery blades, round handle. Light, rotatable and handy for precise cutting



**ERGOPLANT  
DX305R**

150 mm, 6"

Microblade holder with round handle, for precision incisions



**BB367R**

Microsurgery blade sterile packed  
Sales unit:  
**PAK** = Package of 10 pieces



**BB369R**

Microsurgical blade stainless steel, sterile packed  
Sales unit: **ST** = package of 1 piece



**BB063R**

145 mm, 5<sup>3</sup>/<sub>4</sub>"

Scalpel handle for sterile scalpel blades



**BB511**

Fig. 11

**BB515**

Fig. 15

**BB515-C**

Fig. 15-C


Sterile scalpel blade, vacuum heat treated, sterile packed  
Sales unit: **PAK** = package of 100 pieces



**BM003R**

180 mm, 7"

DUROGRIP® Needle holder for suture material 4/0 - 6/0

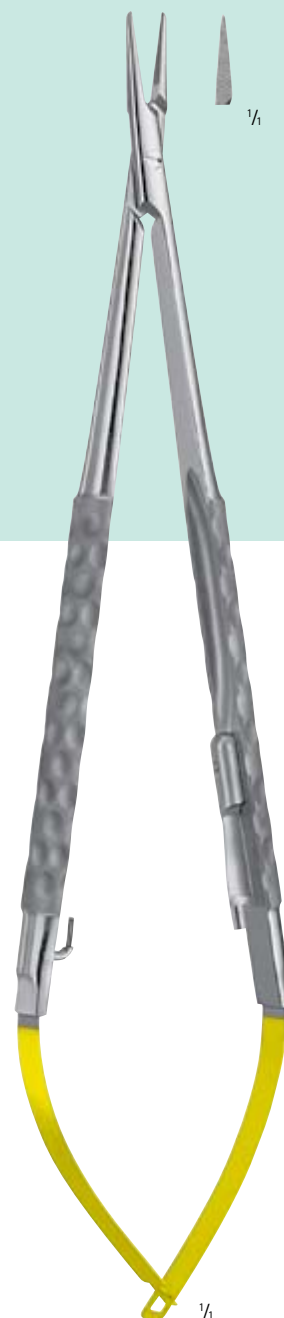


**FD258R**

180 mm, 7"

DIADUST® needle holder with carbide inserts for suture material 6/0 - 8/0

# for periodontal recessions

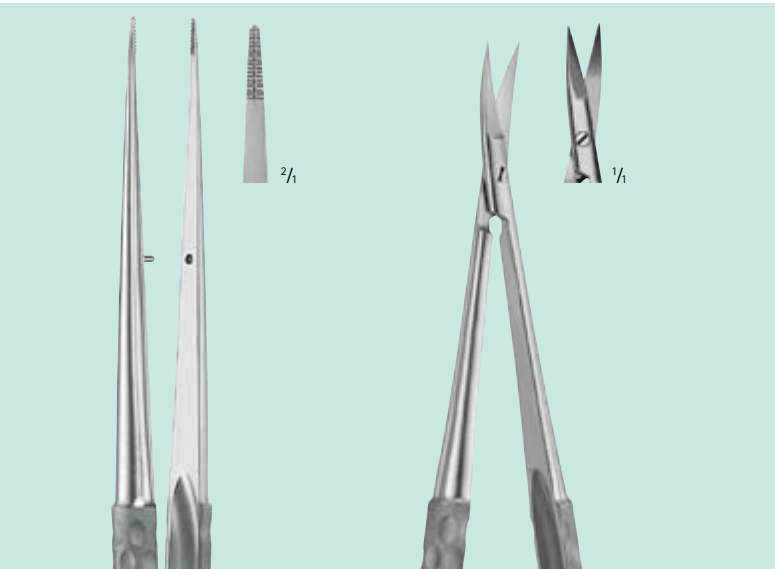


**COOLEY  
DX303R**  
180 mm, 7"  
Atraumatic  
micro forceps

**FD252R**  
180 mm, 7"  
Micro ligature scissors, curved,  
one blade serrated



**JG381R**  
274 x 172 x 30 mm  
Tray, including 3 x JG391  
silicone bars, suitable for  
dental container JN092



Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany  
Phone +49 7461 95-0 | Fax +49 7461 95-26 00 | [www.aesculap-dental.com](http://www.aesculap-dental.com)

Aesculap – a B. Braun company

The main product trademark 'Aesculap' and the product trademarks 'Sterilit', 'DIADUST' and 'DUROGRIP' are registered trademarks of Aesculap AG.

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products. Reproduction in any form partial or otherwise is not permitted.

Brochure No. D12102

0214/0.25/2