

CE Technical Documentation
SMS Surgical Gown Products
 -- Product Introduction

Documentation No.	CE14-01A
Version	A/0
Effective Date	2016/07/27
Page	1/3

1. Products Information

1.1 Products category

Products	UMDNS Code	Classification Rules	Classification
SMS Surgical Gown	11901	Rule 1	I


1.2 Product description (Raw material, Component and Structure)

The SMS surgical gown products consists of several components, including hydrophobic SMS material, a hook and loop, white knitted material, yellow pure polyester.

All gowns are easy to don and include a hook and loop to close the neck and ties for closing the back of the gown. They include functional cuffs constructed of a white knitted material.

All components are assembled together into a whole surgical gown.

1.3 Product picture and structure

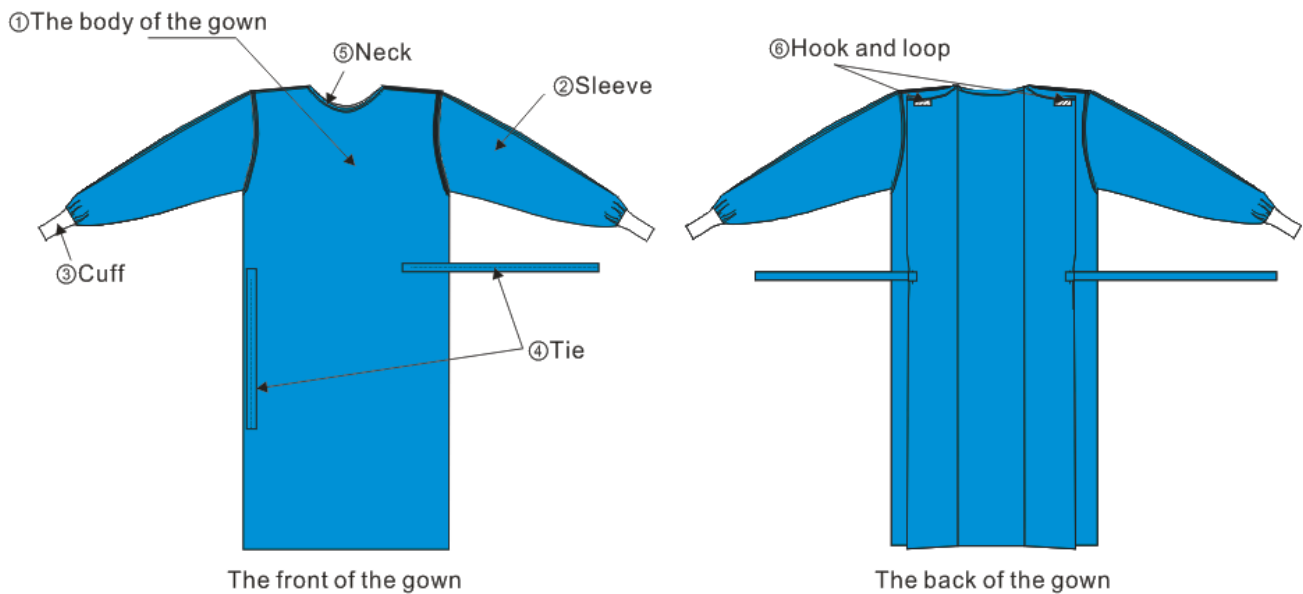

CE Technical Documentation
SMS Surgical Gown Products
 -- Product Introduction

Documentation No.	CE14-01A
Version	A/0
Effective Date	2016/07/27
Page	2/3

Component Name	Material
①The body of the gown	Hydrophobic SMS material 43g, with alcohol-repellent
②Sleeve	Hydrophobic SMS material 43g, with alcohol-repellent
③Cuff	White pure polyester
④Tie	Hydrophobic SMS material 43g, with alcohol-repellent
⑤Neck	Pure polyester
⑥Loop and hook	Mixing polyester with nylon

1.4 Basic size

No.	Product	Size
1	SMS Surgical Gown 43g, with alcohol-repellent	L 126x155cm
2		XL 131x160cm
3		XLL 151x160cm
4		XXL 151x170cm

2. Intended Purpose/Use

The standard disposable SMS surgical gowns are items of surgical apparel intended to be worn by health care professionals to help protect both the patient and health care worker from transfer of microorganisms, body fluids and particulate matter.

All SMS surgical gowns are generously cut in the chest and sleeves for greater comfort and freedom of movement. Each SMS surgical gown features a very long hook & loop neck closure for an easily adjustable, gap-free custom fit. Wide range of sizes and lengths satisfies the needs of every clinician, with high-tech breathable fabrics keeping the wearer comfortable and dry, even during long procedures.